

Pipette filler, Bibbette

Azlon

- This pipette filler offers incredible control.
- The speed of operation is similar to electronic units many times its price.
- A novel PTFE and Silicone rubber valve assembly controls the vacuum
- By moving the thumb lever, pipettes from 0.1 ml up to 100 ml can be easily filled
- Blow out pipettes are catered for by the depression of a small thumb operated rubber bulb.
- A replaceable PTFE membrane filter ensures liquid is not accidentally drawn into the unit (come in packs of 10).
- The filter and pipette holding cone can be sterilised by autoclaving.
- Each Bibbette comes with full instructions and two membrane filters

Cat. No.	Code	Description
PF/13002	PFL001N	Bibbette pipette controller, Plain
PF/13004	PFL001B	Bibbette pipette controller, Blue
PF/13014	PFL002	Spare membrane filters pack of 10
PF/13016	PFL004	Spare valve assembly



PF/130

Pipette fillers, rubber

Moulded rubber bulb approx. 60ml capacity, with three glass ball valves which can be controlled by finger pressure. No metal parts to corrode and can be used for all liquids except those which attack rubber.

Cat. No.	Colour
PF/96112	Red
PF/96114	Black
PF/96116	Yellow
PF/96118	Blue




PF/961

Pipette fillers

- Three sizes, 2ml, 10ml and 25ml
- Easy to use design
- Elastic step tip and drain valve

Cat. No	Code	Capacity	Colour
PF/13118	R622	Up to 2ml	Blue
PF/13120	R623	Up to 10ml	Green
PF/13122	R624	Up to 25ml	Red

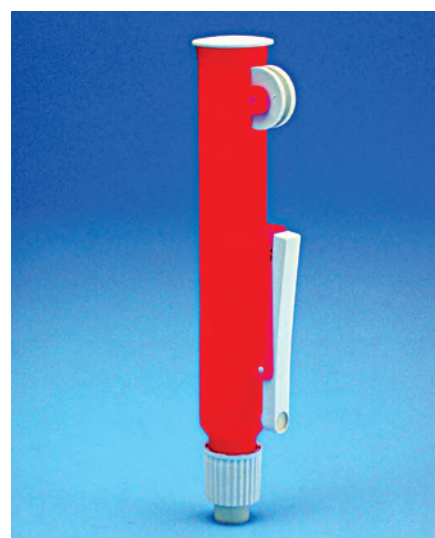

PF/131

Pipette Controller, Bibbyjet Pro, rechargeable

Cole-Parmer

The Bibbyjet is powerful, lightweight and ergonomically designed. It fits easily in the hand and is comfortable to use all day.

- LED shows when battery is running low
- The Bibbyjet has a venting feature. The air sucked from the pipette does not go through the pump but is discharged out of the front of the instrument. This avoids corrosion of the internal mechanism by the damp air and vents any fumes away from the operator.
- Mode selector switch for gravity delivery or blow out (with motor power)

A filter and safety valve protect against contamination of the unit and the pipette holder is fully autoclavable. The new Bibbyjet is powerful enough to fill a 50ml pipette in just 10 seconds and the batteries have the power to give a full day's work without recharging. Supplied complete with 2 replacement membrane filters and wall holder.

Cat. No.	Code	Description
PC/15002	PC2000	Bibbyjet Pro, complete with UK charger
Accessories		
PC/15004	PC1011	Silicone adaptor
PC/15006	PC1012	0.2µm filter, sterile, pack of one
PC/15008	PC1013	0.2µm filter, non-sterile, pack of 10
PC/15010	PC1014	3µm filter, non-sterile, pack of one
PC/15012	PC1015	3µm filter, non-sterile, pack of 10
PC/15014	PC1016	Wall support



PC/15002