



FB250/3SC

## Filter flasks, Buchner filter flasks

Features a unique 'push in' connector which allows tubing to be fitted easily, quickly and in complete safety. Tubing is attached to the plastic connector which is then simply inserted into a socket on the side of the flask and tightened to give a firm vacuum seal. This system offers major benefits over traditional glass side arm connectors as it is:

- Safe: Less risk of breakage and personal injury.
- Easy to use: Simply push in connector and tighten screw cap
- Effective: Provides a firm leak-free seal.
- Economical: Risk of breakage and cost of replacement greatly reduced.
- Durable: Connector components are chemically resistant and are autoclavable.

Cat. No.	Nominal cap, ml	Socket size	Overall height, mm	Approx. dia., mm
FB100/2SC	100	19/26	105	70
FB100/3SC	100	24/29	105	70
FB250/3SC	250	24/29	140	85
FB250/4SC	250	29/32	145	85
FB500/3SC	500	24/29	172	105
FB500/4SC	500	29/32	172	105
FB1L/3SC	1000	24/29	220	135
FB1L/4SC	1000	29/32	220	135




FC7/23

## Fractionating columns, plain

Borosilicate glass with ground-glass joints. Designed for use with various types of column packing e.g. Raschig rings.

Cat. No.	Approx. overall length, mm	Approx. surface area x 10 <sup>-2</sup> m <sup>2</sup>	Approx. effective length, mm	Approx. overall height, mm	Cone size	Socket size
FC7/23	365	2	250	365	24/29	24/29
FC7/43	615	4	500	615	24/29	24/29

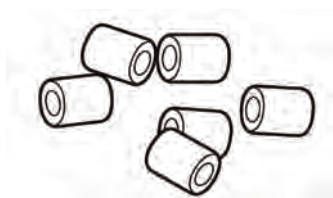

FC8/09

## Column packing

Packing, for use with above plain fractionating columns.

Manufactured from borosilicate glass for contamination-free performance.

Available as Raschig rings in 1 litre bags.

Raschig rings; Short lengths of glass tubing.

Cat. No.	Approx. dia, mm	Approx. length, mm	Description
FC8/09	9	9	Raschig rings



EX13/26

## Packing supports

Packing support, for Raschig rings within plain fractionating columns.

Support is passed through socket of column and rests at the constriction above the cone.

Manufactured from chemically resistant borosilicate glass.

Cat. No.	Diameter, mm	To t socket size
EX13/12	16	19/26
EX13/23	21	24/29
EX13/25	29	34/35
EX13/26	35	40/38